

# ClimaZone

- Lightweight cushion performs at bariatric levels up to 254 kg
- Temperature regulating with CoolOver cover
- For moderate skin protection and mild to moderate postural support needs, contoured version available for high to maximum postural support needs


Lightweight cushion that allows for air circulation to reduce sweat and moisture. Revolutionary construction that reduces pressure using air suspended fibres that promote airflow and integrated CoolOver fabric to regulate skin temperature. It's the go-to-cushion to manage heat and moisture at the seated surface - also offered for bariatric size and weight capacities.

## Cover:

- 1. CoolOver Cover** to regulate skin temperature to prevent skin breakdown, drives away excess heat from the skin, facilitates heat retention when needed, maintains optimal skin temperature
- 2. Microclimatic Cover:** 2-way stretch spacer fabric to promote circulation and reduces skin insulation



## Foam:

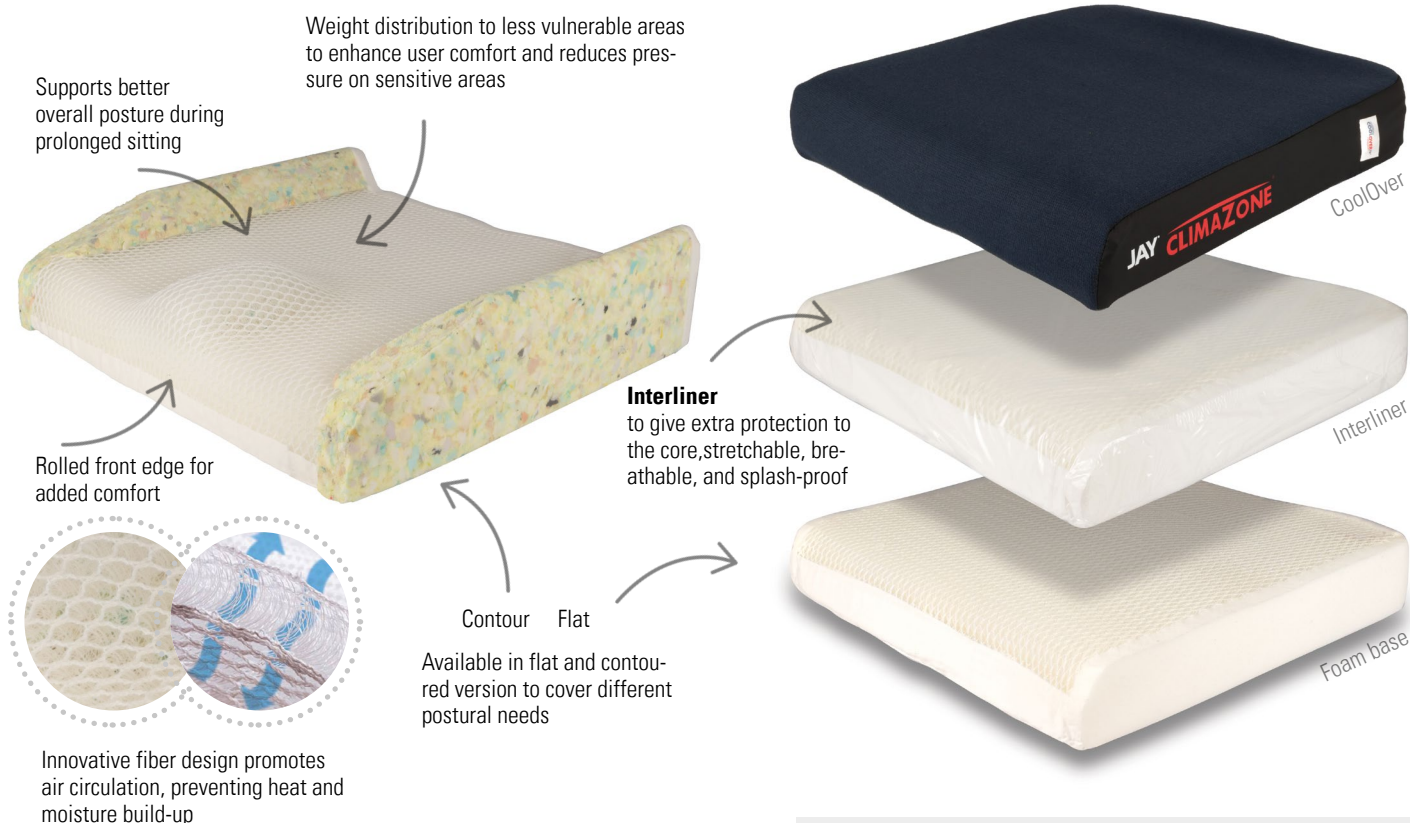

Scan the QR code

QR codes require adaptation



website




order form

These sizes are available, please see your order form for modifications or special sizes.

|                  |                 |
|------------------|-----------------|
| Weight cushion   | 1,00 kg         |
| Max. user weight | 254 kg          |
| Height           | 50 mm and 80 mm |
| Width            | 300 - 660 mm    |
| Depth            | 300 - 600 mm    |

# JAY Cushion range

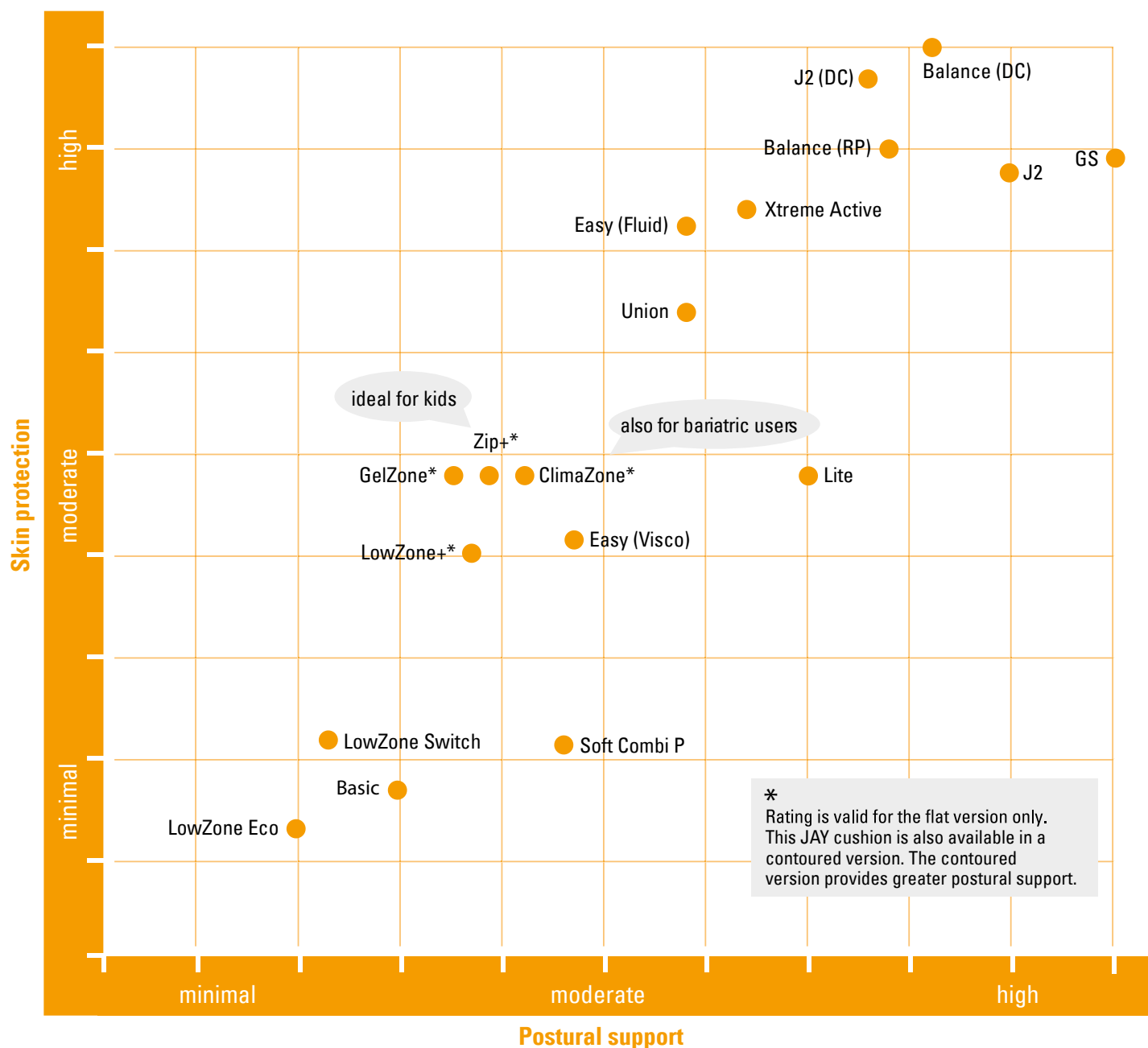

The demands on a seat cushion are as different as people are individual and compromises in clinical use are not acceptable. For this reason, our JAY portfolio offers tailored seating solutions spanning the broad continuum of clinical needs, while maintaining clinical efficacy. All of our JAY cushions combine skin protection with effective postural support and ease of use to the highest degree.

#### DISCLAIMER: FOR PROFESSIONAL USE ONLY.

Sunrise Medical strongly recommends that a clinician, such as a doctor or therapist experienced in seating and positioning, be consulted to determine which cushion and/or back is appropriate for the users' specific needs to prevent injury. Cushions, backs, and positioning supports should only be installed by an authorized Sunrise Medical dealer.

Sunrise Medical GmbH

Kahlbachring 2-4

69254 Malsch/HD

Tel.: +49 (0) 72 53 / 9 80-0

Fax: +49 (0) 72 53 / 9 80-222

E-Mail: [international@sunrise-medical.de](mailto:international@sunrise-medical.de)

© 03.2024



[www.SunriseMedical.eu](http://www.SunriseMedical.eu)